


## TROPHIC CONTINENCE STIMULATOR

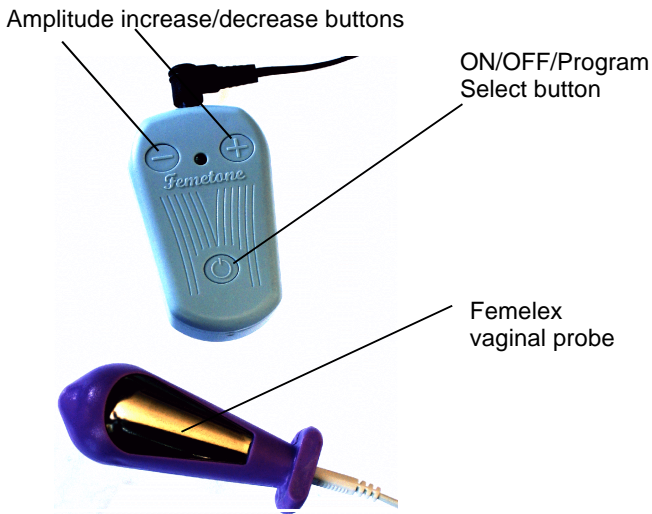

The **Femetone** trophic pelvic floor stimulator is an advanced development in a small, lightweight neuromuscular electrical stimulator (E-stim) for clinic and home use. The programs mimic the natural nerve-muscle communication associated with bladder control. With its simple push-button controls, it is especially easy to use at home.

Unlike other dedicated electrical stimulators, this functional electrical stimulator incorporates automatic sequential modes; this enables the therapist, continence adviser or doctor to more effectively treat and target specific symptoms.

After selecting one of the AUTOMATIC PROGRAM MODES, the muscle is activated through a series of different levels of metabolic activity, from aerobic to voluntary assisted maximum contractions, ending the cycle with a cooling off phase.

### **Automatic Sequential Frequency Programs**

1. Strengthening weak muscles..... 40 minutes
2. Active muscle exercise..... 25 minutes
3. Activating weak atrophied muscles..... 60 minutes

Conditions of stress and urge incontinence often coexist with bladder instability, overactive bladder (OAB) syndrome, sexual dysfunction and pelvic pain to varying degrees. There are many advantages from using a trophic pelvic floor stimulator. **Femetone** produces multiple reflex effects, including restoring urinary control and treating an overactive detrusor. Improvements in urethral closure are a progression towards

### RECOMMENDED FOR TREATMENT OF

- Genuine stress urinary incontinence
- Stress / Urge mixed incontinence
- Unstable bladder
- Symptom-free pelvic floor weakness
- Loss of normal vaginal sensation
- First degree prolapse

### FEATURES

- Selectable programs for muscle tone and muscle strengthening
- Sequential frequency muscle activation

### Femetone Specifications

- Single channel output
- Asymmetrical bi-phasic net zero DC
- 80 milli-amps amplitude, max
- Multi mode operation
- Automatic sequential frequencies
- On cycle 4 seconds
- Off cycle 4 seconds
- Pulse width 220 F second
- Frequencies 3 to 40 pulses per second
- 3-volt CR2032 lithium battery

### Kit includes:

- FT-80 Femetone stimulator
- S31 ISIS Femelex vaginal probe
- S97 ISIS Electrode lead wire
- Instruction manual - English, German, French
- Extra Batteries