


Nuromed

Femetone

TROPHIC CONTINENCE STIMULATOR

Amplitude increase/decrease buttons

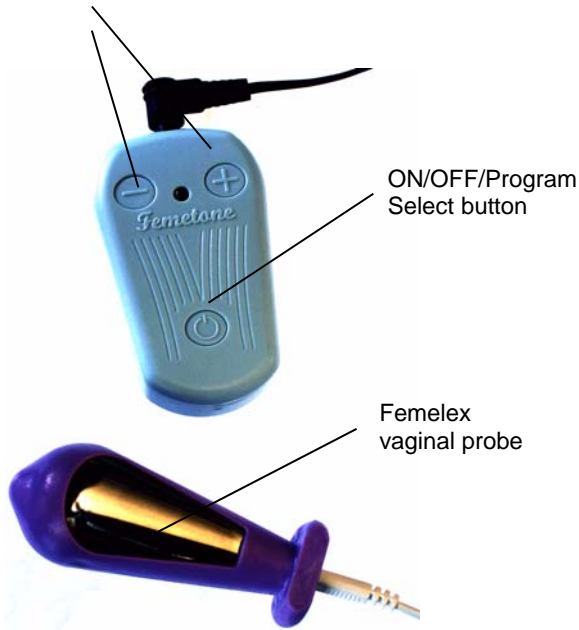

RECOMMENDED FOR TREATMENT OF

- Genuine stress incontinence
- Stress / Urge mixed incontinence
- Unstable bladder
- Symptom free pelvic floor weakness
- Loss of normal vaginal sensation
- First degree prolapse

FEATURES

- Selectable programs for muscle tone and muscle strengthening
- Sequential frequency muscle activation
- Touch button controls
- Both vaginal and anal electrodes available

The *Femetone* trophic continence stimulator is an advanced development in a small, lightweight neuromuscular electrical stimulator. The programs mimic the natural nerve-muscle communication associated with continence control. With its simple push-button controls, it is especially easy to use at home.

Unlike other dedicated electrical stimulators it incorporates automatic sequential modes; this enables the therapist, continence adviser or doctor to more effectively treat and target specific symptoms.

After selecting one of the **AUTOMATIC PROGRAM MODES**, the muscle is activated through a series of different levels of metabolic activity, from aerobic to voluntary assisted maximum contractions, ending the cycle with a cooling off phase.

Automatic Sequential Frequency Programs

- | | |
|--------------------------------------|------------|
| 1. Strengthening weak muscles | 40 minutes |
| 2. Active muscle exercise | 25 minutes |
| 3. Activating weak atrophied muscles | 60 minutes |

Conditions of stress and urge a incontinence often coexist with bladder instability and pelvic pain to varying degrees. There are many advantages of using trophic stimulator. *Femetone* produces multiple reflex effects, including restoring urinary control and treating an overactive detrusor. Improvements in urethral closure are a progression towards complete continence.

Femetone Specifications

- Single channel
- Asymmetrical bi-phasic net zero DC
- 80 milli-amps amplitude, max
- Multi mode operation
- Automatic sequential frequencies
- On cycle 4 seconds
- Off cycle 4 seconds
- Pulse width 200 F second
- Frequencies 5 to 40 pulses per second
- 3 volt CR2032 lithium battery

Kit includes:

- Femetone stimulator
- S31 ISIS Femelex vaginal probe
- S97 ISIS Electrode lead wire
- Instruction manual

BIOMATION

335 Perth Street
P.O. Box 156
Almonte, Ontario K0A 1A0

Tel: (613) 256-2821
Toll Free: 1-888-667-2324
Fax: (613) 256-5872
E-Mail: info@biomation.com
Web: www.biomation.com/inco